

EMERGENCY VACCINE TRANSPORT

DURING AN EMERGENCY WHERE VACCINES HAVE TO BE TRANSFERRED TO ANOTHER LOCATION, PLEASE FOLLOW THESE STEPS:

- NOTIFY R.E. OF TRANSFER AND REASON
- CREATE AN ITEMIZED LIST OF VACCINES BEING TRANSFERRED, THE LIST SHOULD READ AS FOLLOWS: VACCINE NAME + PEDI OR ADULT, QUANTITY LOT# AND EXPIRATION DATE.
- WHEN TRANSFER HAS BEEN COMPLETED EMAIL VACCINE LIST TO R.E., AND THE PIN # WHERE VACCINES WERE TRANSFERRED TO, ALONG WITH THE VSTR REPORT IF ANY.
- MAINTAIN A TEMPERATURE LOG, DOCUMENT TEMPERATURES EVERY 10 MINUTES IF VACCINES ARE NOT IMMEDIATELY PLACED IN THE BACK UP UNIT TRANSFERRED TO.

EMERGENCY VACCINE TRANSPORT - PLEASE CONTACT R.E. IMMEDIATELY IF THIS OCCURS

REQUIRED MATERIALS FOR VFC/ASN VACCINE TRANSPORT



COOLER/S LARGE ENOUGH TO ACCOMMODATE CLINIC'S TYPICAL SUPPLY OF REFRIGERATED AND FROZEN VACCINES.
CAN USE BOXES FROM MANUFACTURER'S IF AVAILABLE
DO NOT USE SOFT - SIDED COLLAPSIBLE COOLERS.



FOR INSULATION
2 - CORRUGATED CARDBOARDS CUT TO FIT COOLER'S INTERIOR



FOR INSULATION
2 - LAYERS OF 1" THICK BUBBLE WRAP CUT TO FIT COOLER'S INTERIOR





**PLACE A COPY OF YOUR EMERGENCY PLAN EITHER INSIDE OR ON TOP OF COOLER.
- MAKE SURE PLAN IS UPDATED YEARLY FROM DATE OF LAST REVIEW.**



MAKE SURE THAT YOU HAVE ENOUGH FROZEN BOTTLES TO LINE BOTH TOP AND BOTTOM OF COOLERS THAT WILL BE USED DURING TRANSPORT.



CONDITION FROZEN WATER BOTTLES - PLACING BOTTLES IN SINK FILLED WITH SEVERAL INCHES OF WARM WATER OR UNDER RUNNING TAP -

BOTTLES ARE CONDITIONED WHEN THE ICE EASILY SPINS INSIDE OF BOTTLE - WIPE BOTTLES OFF AND PLACE IN COOLER.



PACKING COOLER FOR EMERGENCY TRANSPORT OF REFRIGERATED VACCINES

LINE
CONDITIONED
BOTTLES ON
BOTTOM OF
COOLER



1. Conditioned Water Bottles

PLACE
CARDBOARD
ON TOP OF
BOTTLES



2. Cardboard Sheet

ADD LAYER OF
BUBBLE WRAP
ON TOP OF
CARDBOARD




3. Bubble wrap, packing foam, or Styrofoam™

PLACE TEMPERATURE READY PROBE
IN CENTER ENSURING THE CORD
REACHES TO OUTSIDE OF COOLER TO
PLUG IN LOGGER, ADD VACCINES.



4. Vaccines, Diluents, and Temperature Monitoring Device Probe

AFTER VACCINES HAVE BEEN ADDED
FOLLOW STEPS 5 - 7

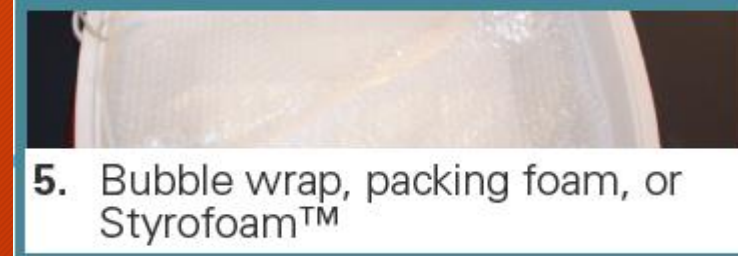

5. Bubble wrap, packing foam, or Styrofoam™

ADD BUBBLE WRAP ON TOP OF
VACCINES



6. Cardboard Sheet

LAY CARDBOARD ON TOP OF BUBBLE
WRAP



7. Conditioned Water Bottles

ADD THE REST OF THE PREPPED BOTTLES AND
CLOSE LID. PLUG IN LOGGER AND POWER UP;
MAKE SURE LOGGER IS CONFIGURED FOR
REFRIGERATED TEMPERATURES.

TANSFERRING FROZEN VACCINES

GREENBOX IS AN IDEAL WAY TO TRANSFER FROZEN VACCINES. NO PREP-TIME, EASY TO PACK, LIGHTER WEIGHT AND MAINTAINS TEMPERATURE SENSITIVE VACCINES.

PLEASE ENSURE YOU HAVE ON HAND PROPER CONTAINERS FOR TRANSFERRING FROZEN VACCINES.

